

Figure 1. Agarose Gel Electrophoresis of DNA isolated from whole blood using MagaCell™ or Qiagen QIAamp DNA Mini Kit, showing high molecular weight DNA isolated by both techniques.

Lane 1: 1 Kb DNA Ladder

Lane 2: Calf thymus DNA Control

Lanes 3, 5, 7, 9, and 11: DNA isolated by MagaCell

Lanes 4, 6, 8, 10, and 12: DNA isolated by QIAamp

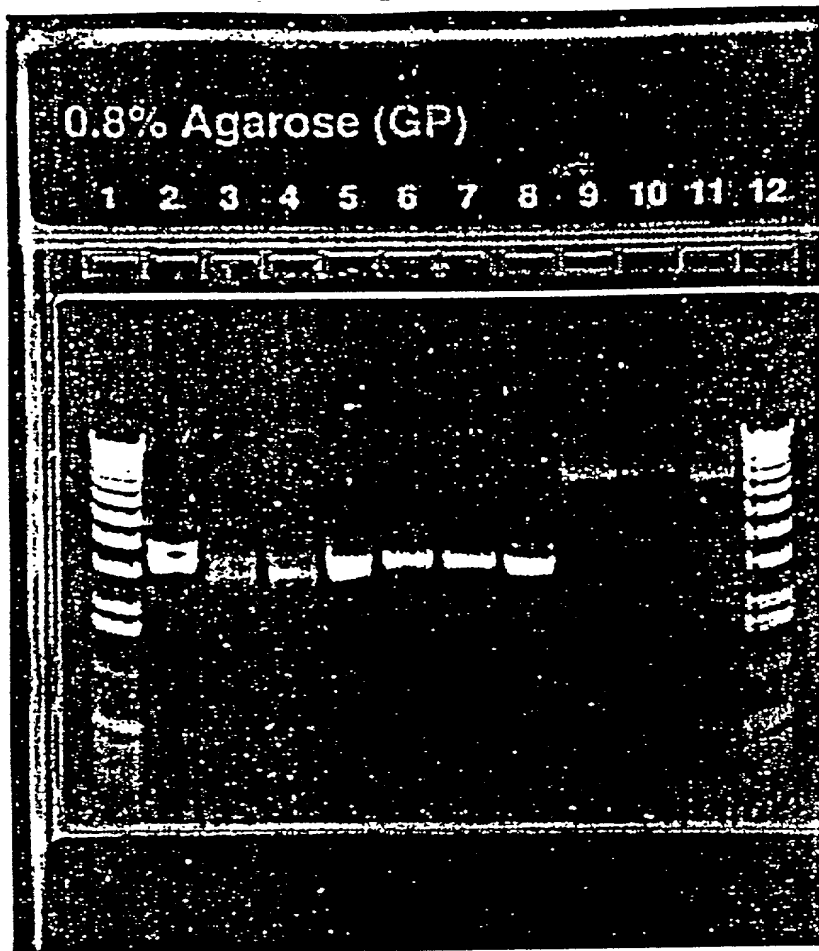


Figure 2. Agarose Gel Electrophoresis of Plasmid DNA isolated from bacterial cell lysates, using MaCell™ or Qiagen QIAprep Miniprep Kit. Two different sizes of high quality plasmid DNA were isolated by both methods.

Lanes 1 and 12: 1 Kb DNA Ladder

Lane 2: Plasmid DNA PBA117 Control

Lanes 3, 4, 6, and 7: Plasmid DNA PBA117 isolated by MagaCell


Lanes 5 and 8: Plasmid DNA PBA117 isolated by QIAprep Miniprep

Lanes 9 and 10: Plasmid DNA PBA8 isolated by MagaCell

Lane 11: Plasmid DNA PBA8 isolated by QIAprep Miniprep

05031 2522660

**Figure 3. Real Time RT-PCR Quantitation of MS2 Viral
RNA Isolated by MagaCell™ or RNeasy Kit**

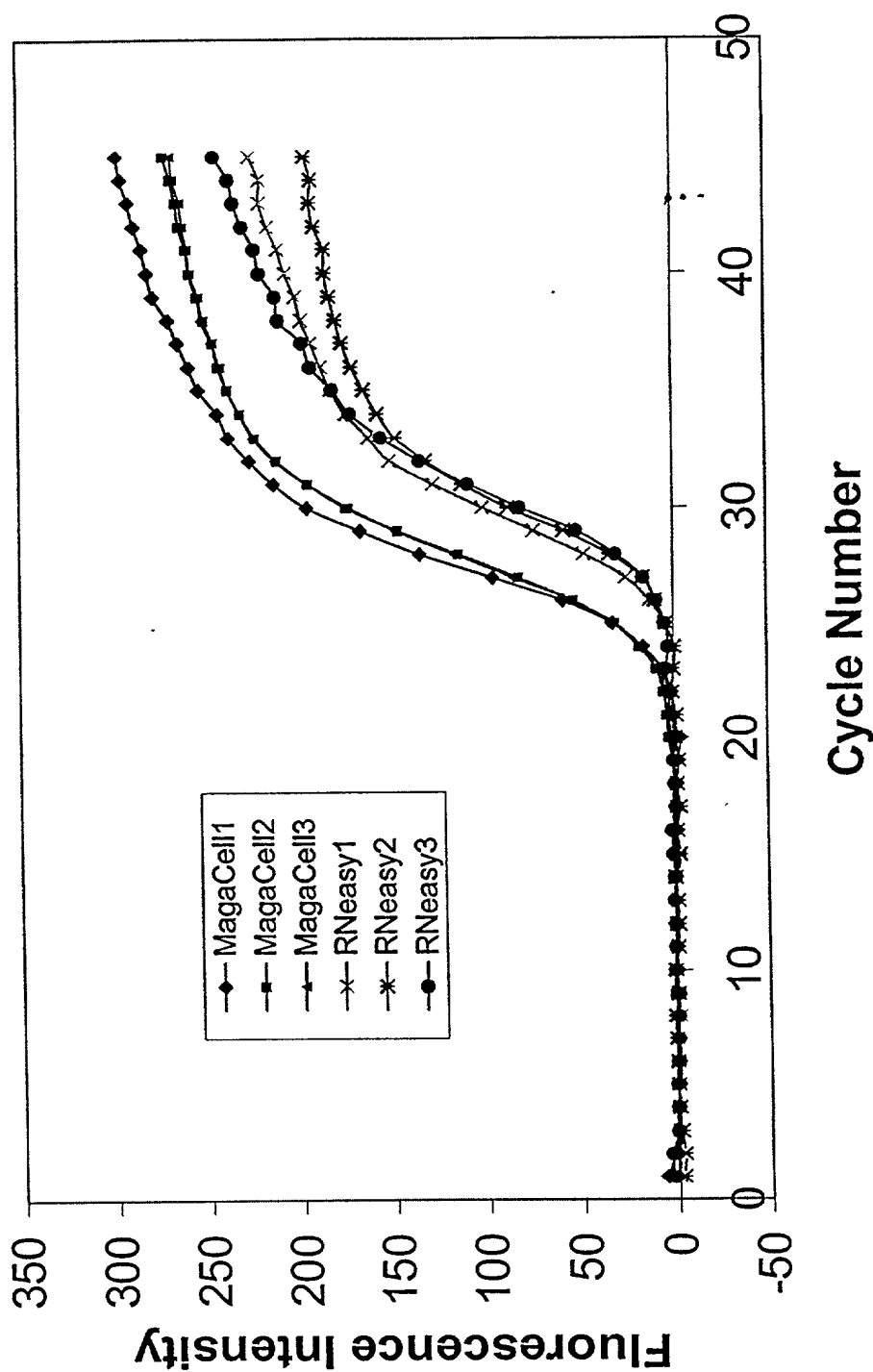

Significant increase in fluorescence intensity was observed after 26 cycles of amplification when MagaCell-purified RNA was used as compared to 29 cycles when RNeasy-purified RNA was used.